PTPN23 Rabbit pAb

Catalog No.: A4595



Basic Information

Observed MW

179kDa

Calculated MW

179kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Mouse, Rat

Background

This gene encodes a member of the non-receptor type protein-tyrosine phosphatase family. The encoded protein may be involved in the regulation of small nuclear ribonucleo protein assembly and pre-mRNA splicing by modifying the survival motor neuron (SMN) complex. The encoded protein additionally plays a role in ciliogenesis and is part of endosomal sorting complex required for transport (ESCRT) pathways. This gene may serve a tumor suppressor function.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

Immunogen Information

Gene IDSwiss Prot
25930
Q9H3S7

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1230-1450 of human PTPN23 (NP_056281.1).

Synonyms

HDPTP; HD-PTP; NEDBASS; PTP-TD14; PTPN23

Contact

<u>a</u>	400-999-6126
\bowtie	cn.market@abclonal.com.cn
\odot	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

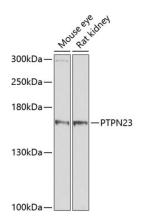


Western blot analysis of various lysates using PTPN23 Rabbit pAb (A4595) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins: $25\mu g$ per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 5s.

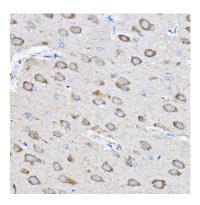

Immunohistochemistry analysis of paraffinembedded Rat brain using PTPN23 Rabbit pAb (A4595) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

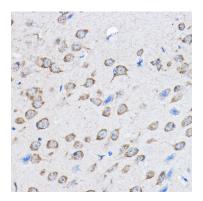

Immunohistochemistry analysis of paraffinembedded Mouse brain using PTPN23 Rabbit pAb (A4595) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.