# ABclonal www.abclonal.com

## Galectin 1/LGALS1 Rabbit mAb

Catalog No.: A4732 Recombinant

## **Basic Information**

#### **Observed MW**

14kDa

## **Calculated MW**

15kDa

## Category

Primary antibody

## **Applications**

ELISA,WB,IHC-P,IF/ICC

#### **Cross-Reactivity**

Human, Mouse, Rat

#### CloneNo number

ARC1110

## **Background**

The galectins are a family of beta-galactoside-binding proteins implicated in modulating cell-cell and cell-matrix interactions. This gene product may act as an autocrine negative growth factor that regulates cell proliferation.

## **Recommended Dilutions**

WB	1:500 - 1:2000	
IHC-P	1:50 - 1:200	
IF/ICC	1:50 - 1:200	

## **Immunogen Information**

Gene ID	Swiss Prot
3956	P09382

#### **Immunogen**

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human Galectin 1/LGALS1 (P09382).

## **Synonyms**

GBP; GAL1; Galectin 1/LGALS1

## Contact

<u>a</u>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\overline{\Box}$	П	www abclonal com cn


## **Product Information**

Source	Isotype	Purification
Rabbit	IgG	Affinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

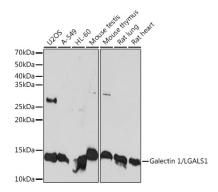


Western blot analysis of various lysates using Galectin 1/LGALS1 Rabbit mAb (A4732) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 3min.

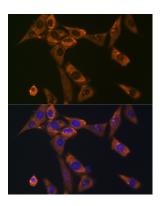

Immunofluorescence analysis of NIH-3T3 cells using Galectin 1/LGALS1 Rabbit mAb (A4732) at dilution of 1:100 (40x lens). Secondary antibody: Cy3 Goat Anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

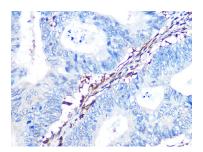

Immunohistochemistry analysis of paraffinembedded Human colon carcinoma using Galectin 1/LGALS1 Rabbit mAb (A4732) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.