Leader in Biomolecular Solutions for Life Science

# **PSMB2** Rabbit pAb

Catalog No.: A5483



### **Basic Information**

Observed MW 23kDa

Calculated MW 23kDa

Category Primary antibody

Applications ELISA,WB,IHC-P,IP

Cross-Reactivity Human, Mouse, Rat

### Background

The proteasome is a multicatalytic proteinase complex with a highly ordered ring-shaped 20S core structure. The core structure is composed of 4 rings of 28 non-identical subunits; 2 rings are composed of 7 alpha subunits and 2 rings are composed of 7 beta subunits. Proteasomes are distributed throughout eukaryotic cells at a high concentration and cleave peptides in an ATP/ubiquitin-dependent process in a non-lysosomal pathway. An essential function of a modified proteasome, the immunoproteasome, is the processing of class I MHC peptides. This gene encodes a member of the proteasome B-type family, also known as the T1B family, that is a 20S core beta subunit. Multiple alternatively spliced transcript variants encoding distinct isoforms have been found for this gene.

### **Recommended Dilutions**

| WB    | 1:500 - 1:1000   |
|-------|--|
| IHC-P | 1:50 - 1:200   |
| IP    | 0.5µg-4µg antibody for<br>200µg-400µg extracts of<br>whole cells |

### **Immunogen Information**

**Gene ID** 5690

Swiss Prot P49721

#### Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-201 of human PSMB2 (NP\_002785.1).

Synonyms HC7-I; PSMB2

### Contact


| 6         | 400-999-6126              |
|-----------|---------------------------|
| $\bowtie$ | cn.market@abclonal.com.cn |
| €         | www.abclonal.com.cn       |


### **Product Information**

**Source** Rabbit **Isotype** IgG Purification Affinity purification


### Storage

Store at -20°C. Avoid freeze / thaw cycles. Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

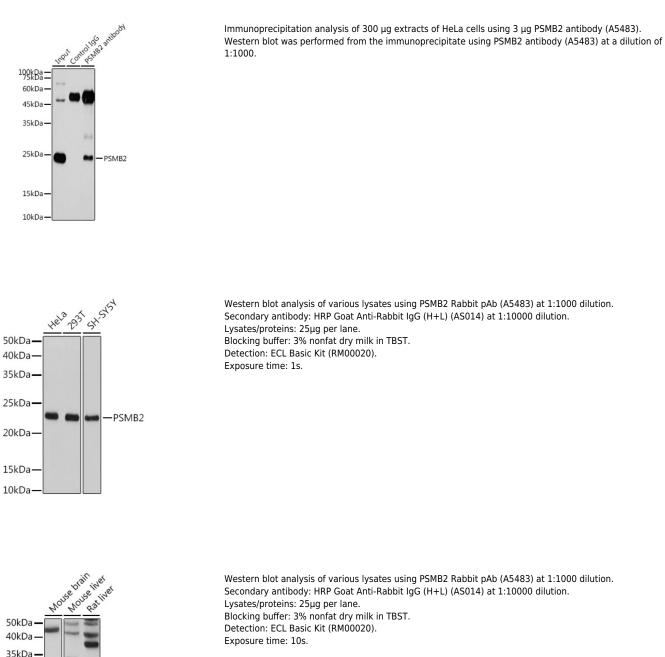

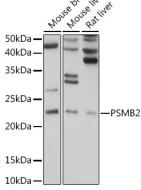

## **Validation Data**

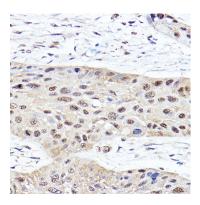

Immunohistochemistry analysis of paraffinembedded Human esophageal cancer using PSMB2 Rabbit pAb (A5483) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.