CLCN1 Rabbit pAb

Catalog No.: A5739



Basic Information

Observed MW

109kDa

Calculated MW

109kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

Background

The CLCN family of voltage-dependent chloride channel genes comprises nine members (CLCN1-7, Ka and Kb) which demonstrate quite diverse functional characteristics while sharing significant sequence homology. The protein encoded by this gene regulates the electric excitability of the skeletal muscle membrane. Mutations in this gene cause two forms of inherited human muscle disorders: recessive generalized myotonia congenita (Becker) and dominant myotonia (Thomsen). Alternative splicing results in multiple transcript variants.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

Immunogen Information

 Gene ID
 Swiss Prot

 1180
 P35523

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 779-988 of human CLCN1 (NP_000074.3).

Synonyms

CLC1; CLCN1

Contact

<u>a</u>	400-999-6126
\bowtie	cn.market@abclonal.com.cn
	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

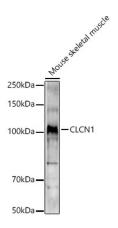

Western blot analysis of lysates from Mouse skeletal muscle, using CLCN1 Rabbit pAb (A5739) at 1:1000 dilution

Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).

Exposure time: 180s.



Immunohistochemistry analysis of paraffinembedded Mouse skeletal muscle using CLCN1 Rabbit pAb (A5739) at dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

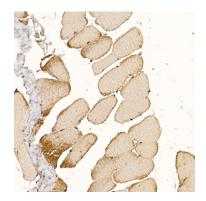

Immunohistochemistry analysis of paraffinembedded Rat skeletal muscle using CLCN1 Rabbit pAb (A5739) at dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.