LALBA Rabbit pAb

Catalog No.: A6233



Basic Information

Observed MW

33kDa

Calculated MW

16kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Mouse, Rat

Background

This gene encodes alpha-lactalbumin, a principal protein of milk. Alpha-lactalbumin forms the regulatory subunit of the lactose synthase (LS) heterodimer and beta 1,4-galactosyltransferase (beta4Gal-T1) forms the catalytic component. Together, these proteins enable LS to produce lactose by transfering galactose moieties to glucose. As a monomer, alpha-lactalbumin strongly binds calcium and zinc ions and may possess bactericidal or antitumor activity. A folding variant of alpha-lactalbumin, called HAMLET, likely induces apoptosis in tumor and immature cells.

Recommended Dilutions

WB 1:500 - 1:2000

IHC-P 1:50 - 1:200

Immunogen Information

Gene ID3906

Swiss Prot
P00709

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 20-142 of human LALBA (NP_002280.1).

Synonyms

LYZG; HAMLET; LALBA

Contact

a		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	Ī	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

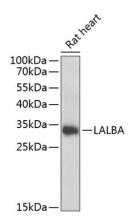


Western blot analysis of lysates from rat heart, using LALBA Rabbit pAb (A6233) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG(H+L)(AS014) at 1:10000 dilution.

Lysates/proteins: 25µg per lane.


Blocking buffer: 3% nonfat dry milk in TBST.

Detection: ECL Basic Kit (RM00020).

Exposure time: 60s.

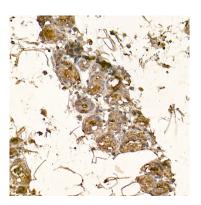

Immunohistochemistry analysis of paraffinembedded Mouse breast using LALBA Rabbit pAb (A6233) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

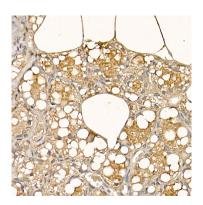

Immunohistochemistry analysis of paraffinembedded Rat breast using LALBA Rabbit pAb (A6233) at dilution of 1:100 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.