# **BCAP29 Rabbit pAb**

Catalog No.: A6335



## **Basic Information**

## **Observed MW**

28kDa

## **Calculated MW**

28kDa

## Category

Primary antibody

## **Applications**

ELISA,WB,IHC-P

#### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

Involved in osteoblast differentiation. Located in membrane.

# **Recommended Dilutions**

**WB** 1:500 - 1:2000

**IHC-P** 1:50 - 1:200

# **Immunogen Information**

**Gene ID**Swiss Prot
55973
Q9UHQ4

## **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 125-241 of human BCAP29 (NP\_061332.2).

## **Synonyms**

B29; BAP29; BCAP29

## Contact

<u>a</u>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
$\overline{\triangle}$	ī	www.ahclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

## **Storage**

Store at -20  $^{\circ}$ C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide,50% glycerol,pH7.3.

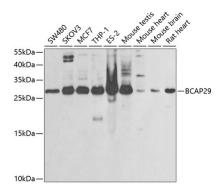


Western blot analysis of various lysates using BCAP29 Rabbit pAb (A6335) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).


Exposure time: 30s.

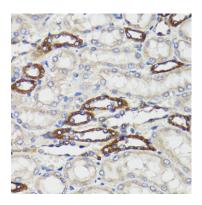

Immunohistochemistry analysis of paraffinembedded Rat kidney using BCAP29 Rabbit pAb (A6335) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

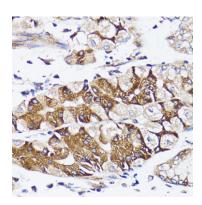

Immunohistochemistry analysis of paraffinembedded Human stomach using BCAP29 Rabbit pAb (A6335) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

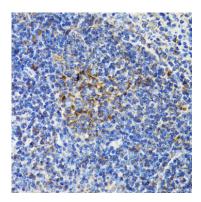

Immunohistochemistry analysis of paraffinembedded Mouse spleen using BCAP29 Rabbit pAb (A6335) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.