EIF2S3 Rabbit pAb

Catalog No.: A6581



Basic Information

Observed MW

51kDa

Calculated MW

51kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

Background

The protein encoded by this gene is the largest subunit of a heterotrimeric GTP-binding protein involved in the recruitment of methionyl-tRNA(i) to the 40 S ribosomal subunit.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:50 - 1:200

Immunogen Information

Gene ID1968

Swiss Prot
P41091

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 1-250 of human EIF2S3 (NP_001406.1).

Synonyms

EIF2; EIF2G; MEHMO; MRXSBRK; eIF-2gA; EIF2gamma; EIF2S3

Contact

<u>a</u>	400-999-6126
\bowtie	cn.market@abclonal.com.cn
\odot	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.01% thimerosal,50% glycerol,pH7.3.

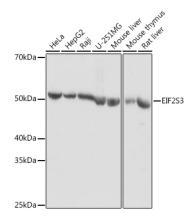


Western blot analysis of various lysates using EIF2S3 Rabbit pAb (A6581) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).


Exposure time: 10s.

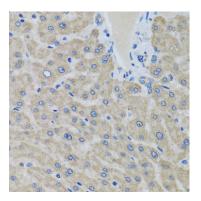

Immunohistochemistry analysis of paraffinembedded Human liver using EIF2S3 Rabbit pAb (A6581) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

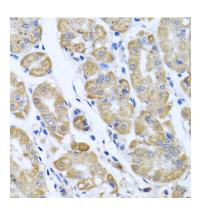

Immunohistochemistry analysis of paraffinembedded Human stomach using EIF2S3 Rabbit pAb (A6581) at dilution of 1:100 (40x lens).Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.

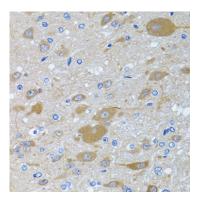

Immunohistochemistry analysis of paraffinembedded Rat brain using EIF2S3 Rabbit pAb (A6581) at dilution of 1:100 (40x lens). Perform microwave antigen retrieval with 10 mM PBS buffer pH 7.2 before commencing with IHC staining protocol.