# **TET3 Rabbit pAb**

Catalog No.: A7612 6 Publications



# **Basic Information**

#### **Observed MW**

179kDa

### **Calculated MW**

194kDa

## Category

Primary antibody

### **Applications**

ELISA, WB, IHC-P, IF/ICC

#### **Cross-Reactivity**

Human, Mouse, Rat

# **Background**

Enables methyl-CpG binding activity and zinc ion binding activity. Involved in histone H3-K4 trimethylation; positive regulation of transcription by RNA polymerase II; and protein O-linked glycosylation. Predicted to be located in cytoplasm and male pronucleus. Predicted to be active in nucleus. Biomarker of esophagus squamous cell carcinoma.

# **Recommended Dilutions**

WB	1:500 - 1:1000
IHC-P	1:50 - 1:200
IF/ICC	1.50 - 1.200

# **Immunogen Information**

Gene ID	Swiss Prot	
200424	043151	

#### **Immunogen**

Recombinant fusion protein containing a sequence corresponding to amino acids 1525-1660 of human TET3 (NP\_001274420.1).

### Synonyms

BEFAHRS; hCG\_40738; TET3

### **Contact**

8		400-999-6126
$\bowtie$	Τ	cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


### **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification

#### Storage

Store at -20°C. Avoid freeze / thaw cycles.


Buffer: PBS with 0.05% proclin300,50% glycerol,pH7.3.

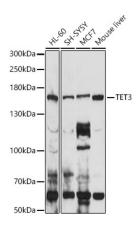


Western blot analysis of various lysates using TET3 Rabbit pAb (A7612) at 1:1000 dilution. Secondary antibody: HRP-conjugated Goat anti-Rabbit lgG (H+L) (AS014) at 1:10000 dilution. Lysates/proteins:  $25\mu g$  per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).

Exposure time: 5s.

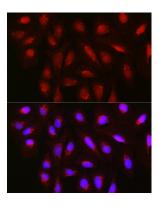

Immunofluorescence analysis of U2OS cells using TET3 Rabbit pAb (A7612) at dilution of 1:100 (40x lens). Secondary antibody: Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007) at 1:500 dilution. Blue: DAPI for nuclear staining.

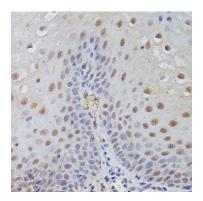

Immunohistochemistry analysis of paraffinembedded Human esophageal using TET3 Rabbit pAb (A7612) at dilution of 1:50 (40x lens). High pressure antigen retrieval performed with 0.01M Citrate Bufferr (pH 6.0) prior to IHC staining.