

# PSD95 Rabbit PolymAb®

Catalog No.: A7889PM

## Basic Information

### Observed MW

95 kDa

### Calculated MW

80 kDa

### Category

Primary antibody

### Applications

WB,IF-F,IHC-P,mIHC,ELISA

### Cross-Reactivity

Mouse, Rat

## Background

This gene encodes a member of the membrane-associated guanylate kinase (MAGUK) family. It heteromultimerizes with another MAGUK protein, DLG2, and is recruited into NMDA receptor and potassium channel clusters. These two MAGUK proteins may interact at postsynaptic sites to form a multimeric scaffold for the clustering of receptors, ion channels, and associated signaling proteins. Multiple transcript variants encoding different isoforms have been found for this gene.

## Recommended Dilutions

<b>WB</b>	1:15000 - 1:65000
<b>IF-F</b>	1:1000 - 1:4000
<b>IHC-P</b>	1:500 - 1:2000
<b>mIHC-P</b>	1:500 - 1:2000
<b>ELISA</b>	Recommended starting concentration is 1 µg/mL. Please optimize the concentration based on your specific assay requirements.

## Contact

☎ | 400-999-6126

✉ | [cn.market@abclonal.com.cn](mailto:cn.market@abclonal.com.cn)

🌐 | [www.abclonal.com.cn](http://www.abclonal.com.cn)

## Immunogen Information

### Gene ID

1742

### Swiss Prot

P78352

### Immunogen

Synthetic peptide. This information is considered to be commercially sensitive.

### Synonyms

MRD62; PSD95; SAP90; SAP-90

## Product Information

### Source

Rabbit

### Isotype

IgG

### Purification


Affinity purification

### Storage




Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS containing 50% glycerol and 0.05% BSA, preserved with proclin300 or sodium azide (as specified on the Certificate of Analysis), pH 7.3.

## Validation Data

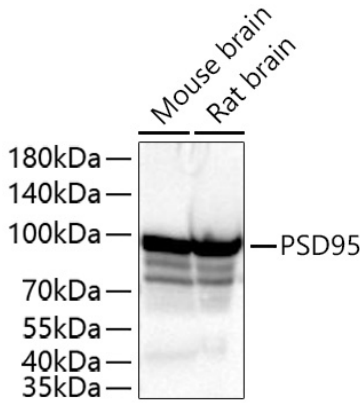




Western blot analysis of various lysates using PSD95 Rabbit PolymAb® (A7889PM) at 1:32000 dilution incubated overnight at 4°C.  
 Secondary antibody: HRP-conjugated Goat anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.  
 Lysates/proteins: 25 µg per lane.  
 Blocking buffer: 3% nonfat dry milk in TBST.  
 Detection: ECL Basic Kit (RM00020).  
 Exposure time: 10 s.

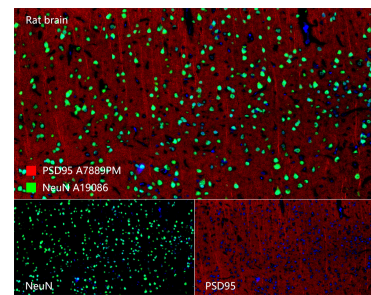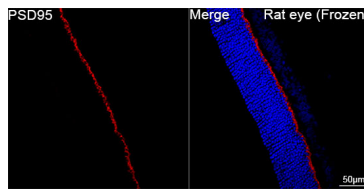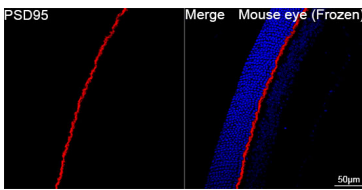

Confocal imaging of frozen sections of Mouse eye tissue using PSD95 Rabbit PolymAb® (A7889PM, dilution 1:4000) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Microwave antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

Confocal imaging of frozen sections of Rat eye tissue using PSD95 Rabbit PolymAb® (A7889PM, dilution 1:4000) followed by a further incubation with Cy3-conjugated Goat anti-Rabbit IgG (H+L) (AS007, dilution 1:500) (Red). DAPI was used for nuclear staining (Blue). Microwave antigen retrieval performed with 0.01M Citrate Buffer (pH 6.0) prior to IF staining. Objective: 40x.

The multiplex IHC analysis on paraffin-embedded Rat brain tissue using the following specific primary antibodies and tyramide signal amplification (TSA) reagents (RK05903) : NeuN Rabbit mAb (A19086, 1:2000) with TSA-TYR-520 (Green), and PSD95 Rabbit PolymAb® (A7889PM, 1:1000) with TSA-TYR-570 (Red). DAPI (Blue) was used for nuclear staining. Prior to multiplex IHC staining, high-pressure antigen retrieval was performed using 0.01M citrate buffer at pH 6.0. The analysis was completed using a 40x objective lens.

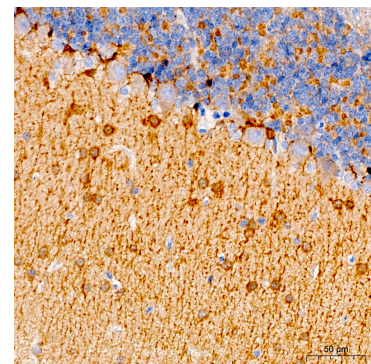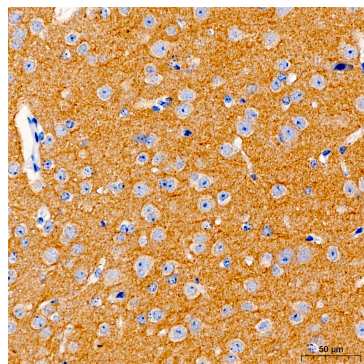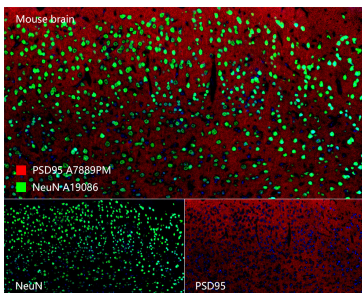

The multiplex IHC analysis on paraffin-embedded Mouse brain tissue using the following specific primary antibodies and tyramide signal amplification (TSA) reagents (RK05903) : NeuN Rabbit mAb (A19086, 1:2000) with TSA-TYR-520 (Green), and PSD95 Rabbit PolymAb® (A7889PM, 1:1000)

Immunohistochemistry analysis of paraffin-embedded Mouse brain tissue using PSD95 Rabbit PolymAb® (A7889PM) at a dilution of 1:800 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

Immunohistochemistry analysis of paraffin-embedded Rat brain tissue using PSD95 Rabbit PolymAb® (A7889PM) at a dilution of 1:800 (40x lens). High pressure antigen retrieval performed with 0.01M Tris-EDTA Buffer (pH 9.0) prior to IHC staining.

## Validation Data

---

with TSA-TYR-570 (Red). DAPI (Blue) was used for nuclear staining. Prior to multiplex IHC staining, high-pressure antigen retrieval was performed using 0.01M citrate buffer at pH 6.0. The analysis was completed using a 40x objective lens.