RBP4 Rabbit mAb

Catalog No.: A8807 Recombinant



Basic Information

Observed MW

22kDa

Calculated MW

21kDa

Category

Primary antibody

Applications

ELISA,WB,IF/ICC

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC1311

Background

This protein belongs to the lipocalin family and is the specific carrier for retinol (vitamin A alcohol) in the blood. It delivers retinol from the liver stores to the peripheral tissues. In plasma, the RBP-retinol complex interacts with transthyretin which prevents its loss by filtration through the kidney glomeruli. A deficiency of vitamin A blocks secretion of the binding protein posttranslationally and results in defective delivery and supply to the epidermal cells.

Recommended Dilutions

WB 1:500 - 1:1000

IF/ICC 1:50 - 1:200

Immunogen Information

Gene ID Swiss Prot 5950 P02753

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 102-201 of human RBP4 (P02753).

Synonyms

RDCCAS; MCOPCB10; RBP4

Contact

<u>a</u>		400-999-6126
\bowtie		cn.market@abclonal.com.cn
\overline{a}	ī	www.ahclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

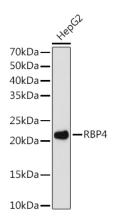


Western blot analysis of extracts of HepG2 cells, using RBP4 antibody (A8807) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Enhanced Kit (RM00021).


Exposure time: 180s.

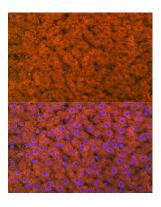

Immunofluorescence analysis of rat liver using RBP4 Rabbit mAb (A8807) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

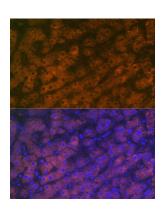

Immunofluorescence analysis of human liver using RBP4 Rabbit mAb (A8807) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

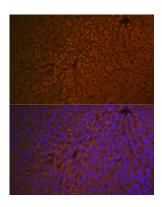

Immunofluorescence analysis of mouse liver using RBP4 Rabbit mAb (A8807) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.