# **CAP1 Rabbit mAb**

Catalog No.: A9330 Recombinant



### **Basic Information**

#### **Observed MW**

46kDa

### **Calculated MW**

52kDa

### Category

Primary antibody

### **Applications**

ELISA,WB,IHC-P

### **Cross-Reactivity**

Human

#### CloneNo number

ARC2752

## **Background**

The protein encoded by this gene is related to the S. cerevisiae CAP protein, which is involved in the cyclic AMP pathway. The human protein is able to interact with other molecules of the same protein, as well as with CAP2 and actin. Alternatively spliced transcript variants encoding different isoforms have been identified.

## **Recommended Dilutions**

**WB** 1:100 - 1:500

**IHC-P** 1:50 - 1:200

## **Immunogen Information**

**Gene ID**Swiss Prot
10487
Q01518

#### **Immunogen**

A synthetic peptide corresponding to a sequence within amino acids 1-100 of human CAP1 (Q01518).

### **Synonyms**

CAP; CAP1-PEN; CAP1

## **Contact**

<b>a</b>		400-999-6126
$\bowtie$		cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


## **Product Information**

SourceIsotypePurificationRabbitIgGAffinity purification


### **Storage**

Store at -20°C. Avoid freeze / thaw cycles.


Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

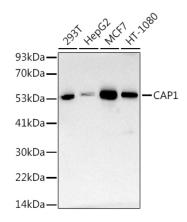

Western blot analysis of various lysates, using CAP1 Rabbit mAb (A9330) at 1:500 dilution.Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.Lysates/proteins: 25ug per lane.Blocking buffer: 3% nonfat dry milk in TBST.Detection: ECL Basic Kit (RM00020).Exposure time: 90s.

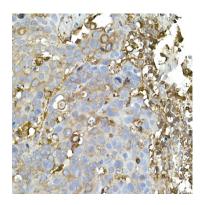

Immunohistochemistry analysis of paraffinembedded Human esophageal cancer using CAP1 Rabbit mAb (A9330) at dilution of 1:50 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

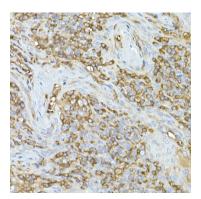

Immunohistochemistry analysis of paraffinembedded Human tonsil using CAP1 Rabbit mAb (A9330) at dilution of 1:50 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.