SAMM50 Rabbit mAb

Catalog No.: A9355 Recombinant



Basic Information

Observed MW

52kDa

Calculated MW

52kDa

Category

Primary antibody

Applications

ELISA,WB,IHC-P

Cross-Reactivity

Human, Mouse, Rat

CloneNo number

ARC2784

Background

This gene encodes a component of the Sorting and Assembly Machinery (SAM) of the mitochondrial outer membrane. The Sam complex functions in the assembly of beta-barrel proteins into the outer mitochondrial membrane.

Recommended Dilutions

WB 1:500 - 1:1000

IHC-P 1:100 - 1:500

Immunogen Information

Gene ID Swiss Prot 25813 Q9Y512

Immunogen

A synthetic peptide corresponding to a sequence within amino acids 300-400 of human SAMM50 (Q9Y512).

Synonyms

OMP85; SAM50; TOB55; TRG-3; CGI-51; YNL026W; SAMM50

Contact

2		400-999-6126
\bowtie		cn.market@abclonal.com.cn
•	T	www.abclonal.com.cn


Product Information

SourceIsotypePurificationRabbitIgGAffinity purification

Storage

Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.02% sodium azide, 0.05% BSA, 50% glycerol, pH7.3.

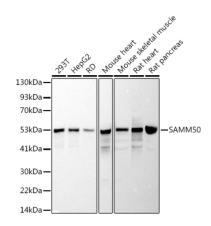


Western blot analysis of various lysates using (A9355) at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.


Lysates/proteins: 25µg per lane.

Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit (RM00020).


Exposure time: 90s.

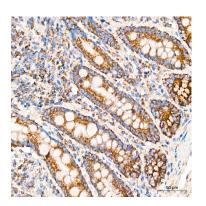


Immunohistochemistry analysis of SAMM50 in paraffin-embedded human colon tissue using SAMM50 Rabbit mAb (A9355) at a dilution of 1:400 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.

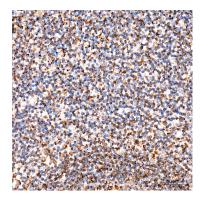


Immunohistochemistry analysis of SAMM50 in paraffin-embedded human tonsil tissue using SAMM50 Rabbit mAb (A9355) at a dilution of 1:400 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.

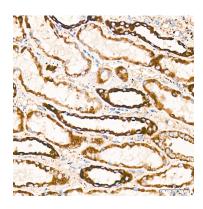


Immunohistochemistry analysis of SAMM50 in paraffin-embedded human kidney tissue using SAMM50 Rabbit mAb (A9355) at a dilution of 1:400 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.

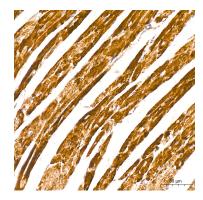

Immunohistochemistry analysis of SAMM50 in paraffin-embedded mouse heart tissue using SAMM50 Rabbit mAb (A9355) at a dilution of 1:400 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.

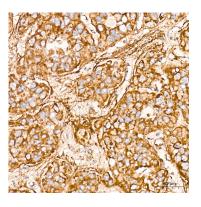

Immunohistochemistry analysis of SAMM50 in paraffin-embedded human liver cancer tissue using SAMM50 Rabbit mAb (A9355) at a dilution of 1:400 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.

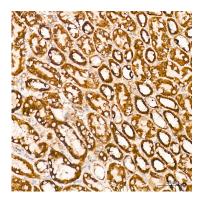

Immunohistochemistry analysis of SAMM50 in paraffin-embedded rat kidney tissue using SAMM50 Rabbit mAb (A9355) at a dilution of 1:400 (40x lens). High pressure antigen retrieval was performed with 0.01 M citrate buffer (pH 6.0) prior to IHC staining.