

β-Tubulin Rabbit pAb

Catalog No.: AC008

100 Publications

Basic Information

Observed MW

55kDa

Calculated MW

50kDa

Category

Loading control antibody

Applications

WB, IHC-P, IF/ICC

Cross-Reactivity

Human, Mouse, Rat

Background

This gene encodes a beta tubulin protein. This protein forms a dimer with alpha tubulin and acts as a structural component of microtubules. Mutations in this gene cause cortical dysplasia, complex, with other brain malformations 6. Alternative splicing results in multiple splice variants. There are multiple pseudogenes for this gene on chromosomes 1, 6, 7, 8, 9, and 13.

Recommended Dilutions

WB	1:500 - 1:1000
IHC-P	1:50 - 1:200
IF/ICC	1:50 - 1:200

Immunogen Information

Gene ID

203068

Swiss Prot

P07437

Immunogen

Recombinant fusion protein containing a sequence corresponding to amino acids 40-260 of human β-Tubulin (NP_006077.2).

Synonyms

M40; TUBB1; TUBB5; CDCBM6; CSCSC1; OK/SW-cl.56

Contact

☎	400-999-6126
✉	cn.market@abclonal.com.cn
🌐	www.abclonal.com.cn

Product Information

Source

Rabbit

Isotype

IgG

Purification


Affinity purification

Storage


Store at -20°C. Avoid freeze / thaw cycles.

Buffer: PBS with 0.05% proclin300, 50% glycerol, pH7.3.


Validation Data

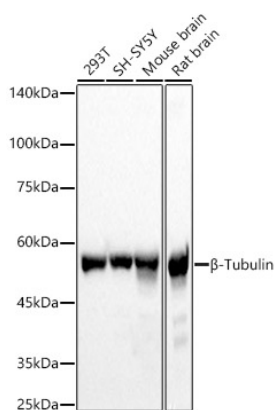


Western blot analysis of various lysates, using β-Tubulin Rabbit pAb (AC008) at 1:2000 dilution.
 Secondary antibody: HRP Goat Anti-Rabbit IgG (H+L) (AS014) at 1:10000 dilution.
 Lysates/proteins: 25μg per lane.
 Blocking buffer: 3% nonfat dry milk in TBST.
 Detection: ECL Basic Kit (RM00020).
 Exposure time: 1s.

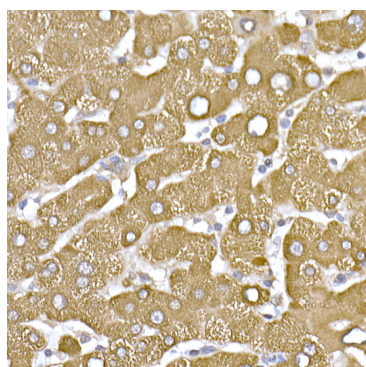


Immunohistochemistry analysis of paraffin-embedded human liver using β-Tubulin Rabbit pAb (AC008) at dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

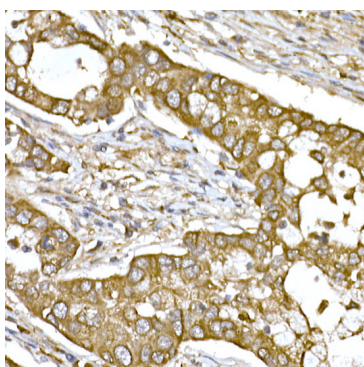


Immunohistochemistry analysis of paraffin-embedded human lung cancer using β-Tubulin Rabbit pAb (AC008) at dilution of 1:200 (40x lens). Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

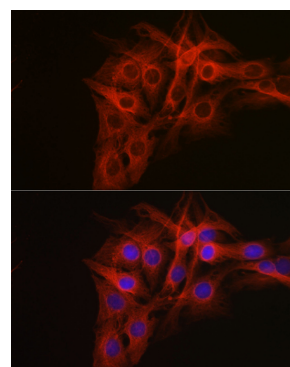

Immunofluorescence analysis of C6 cells using β-Tubulin Rabbit pAb (AC008) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

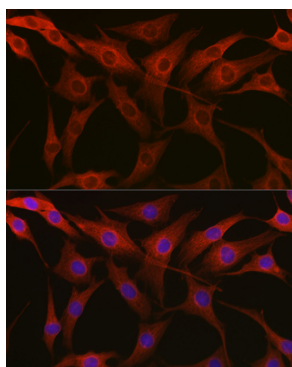

Immunofluorescence analysis of NIH-3T3 cells using β-Tubulin Rabbit pAb (AC008) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.

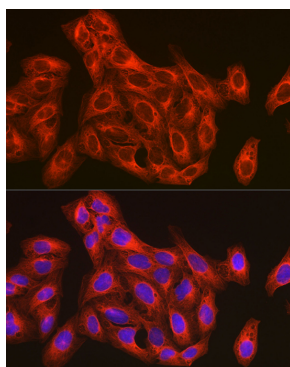

Immunofluorescence analysis of U-2 OS cells using β-Tubulin Rabbit pAb (AC008) at dilution of 1:100 (40x lens). Blue: DAPI for nuclear staining.